

The Accent IP allows users to adjust the heel height up to 5cm with a simple push of a button, giving them complete control over their footwear choices in a low-maintenance package.

A cosmetically appealing foot, the Accent IP is provided with a contoured, soft ankle fairing that locks into the shell and prevents cosmetic buckling when the foot is adjusted at various heights.

This version includes an integrated pyramid for a quick and easy connection to standard 30mm pylon pyramid systems. The Accent's balance of dynamic response and compliance make it a great choice for a vast array of low to moderate intensity activities.

CES CPI F04: Use this code to read the Accent IP clinical evaluation summary on the Steeper website.

Features and Benefits

- User adjustable heel height with up to 5cm range and 13 increments of adjustment
- Easy to operate height adjustment button located in the ankle, available in a choice of 2 locations - medial or lateral
- Ankle fairing provides superior cosmetic finishing
- Suitable for a range of shoe styles




Accommodates various heel heights



Heel height adjustment button Available on Medial or Lateral Side

Cosmetics

The Accent IP's bonded sandal toe foot shell is reinforced with Impact Dispersion Technology™. The Ankle Fairing interlocks with the foot shell and conforms to both Endo and dynamic pylons. The foot shell and ankle fairing are anatomically designed for use with a skin to produce a natural finish.

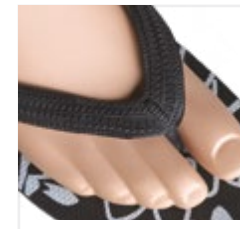

Sandal toe



Ankle fairing



Ordering Information

Trial Period - This foot has a 30 day trial period.

Impact Level Descriptions

Low: Daily activities include mostly level ground walking, moving around the home and the community.

Moderate: Daily activities include up to unlimited walking, climbing stairs and occasional moderate lifting (this does not include running).

Ordering Example: NAC-LCIPL-22N-G2

Part ID	Side	Shell Colour	Mounting	Button Placement	Size	Width**	Firmness Category
NAC	L		IP				G
NAC	R		IP				G

Caucasian	C
Brown	B

Medial	M
Lateral	L

21-28cm

Narrow	N
Wide** 23-26cm only	W

1-3

Firmness Category						
Weight (kg)	0-54	55-63	64-72	73-81	82-100	101-113
Size (cm)	21-28					25-28
Low Impact	1	1	2	2	3	3
Moderate Impact	1	2	2	3	3	3

Spare Parts Ordering Information

Shelltread	Size (cm)
NAC-T ■ N-▲ N	21-28
NAC-T ■ N-▲ W	23-26

- = "L" or "R" for Left or Right
- ▲ = Foot size in cm
- N = Narrow (21-28),
- W = Wide (23-26)





Neutral

Ankle Fairing		
Left	Right	Description
NAC-FLC-2124	NAC-FRC-2124	Caucasian, 21-24cm
NAC-FLC-2526	NAC-FRC-2526	Caucasian, 25-26cm
NAC-FLC-2728	NAC-FRC-2728	Caucasian, 27-28cm
NAC-FLB-2124	NAC-FRB-2124	Brown, 21-24cm
NAC-FLB-2526	NAC-FRB-2526	Brown, 25-26cm
NAC-FLB-2728	NAC-FRB-2728	Brown, 27-28cm

Product Specification

Mounting	Sizes	Weight Limit	Weight of Unit	Build Height	Impact Level	Warranty	Foot Shell
Endo (IP)	21-24cm 25-28cm	100kg 113kg	608g Size 26 inc. foot shell	8.5-9.9cm		2 years	Caucasian or Brown Narrow or Wide width

Technical Information

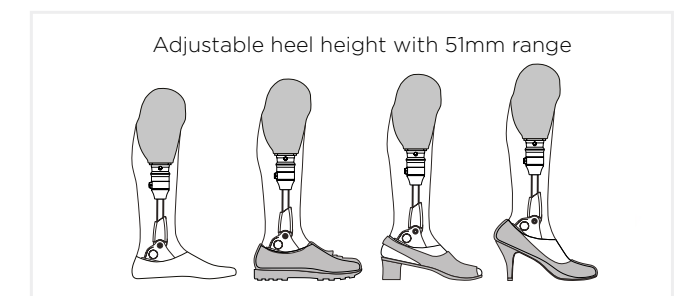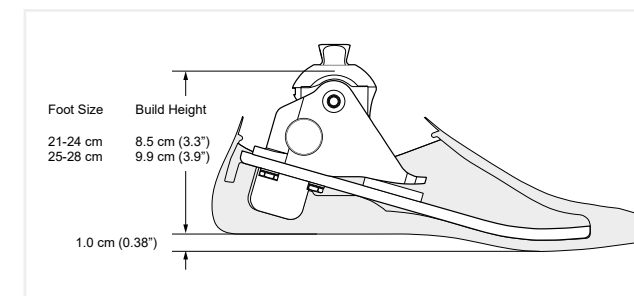

To order contact customerservices@steepergroup.com