

The Odyssey K2 enables the user to walk comfortably and negotiate everyday obstacles such as stairs, curbs and uneven surfaces.

A unique, curved hydraulic cylinder allows independent adjustment of plantar flexion and dorsiflexion via easily accessible dual Stride Control™ dials. The hydraulic ankle, which works in tandem with the Intelliweave® keel, improves the user's posture, knee stability, and comfort, especially when negotiating slopes.


CES CPI F14: Use this code to read the Odyssey K2 clinical evaluation summary on the Steeper website.

Features and Benefits

- Low profile hydraulic foot/ankle combination
- Hydraulic ankle increases ground contact throughout stance, enhancing knee stability while acting to reduce socket interface forces
- Intelliweave® composite foot base provides controlled energy return and a flexible heel
- 12° of smooth hydraulic motion
- Easily accessible, independently adjustable hydraulic valves for tuning in clinic



One pivot angle

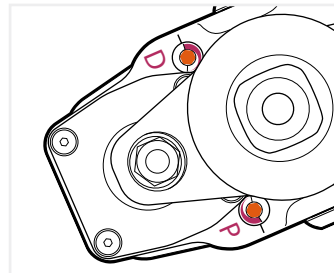

Customisable Stride Control®

Cosmetics



Solid toe foot shell



Ordering Information

Trial Period - This foot has a 30 day trial period.

Impact Level Descriptions

Low: Daily activities include mostly level ground walking, moving around the home and the community.

Ordering Example: NOK2-LCEN-24-G3

Part ID	Side	Shell Colour	Mounting	Size	Firmness Category
NOK2	L		EN		G
NOK2	R		EN		G

Colour	Code
Caucasian	C
Tan	T
Brown	B

Size (cm)	Firmness Category
21-30cm	1-4 See Chart Below

Firmness Category	0-63	64-81	82-100	101-136
Weight (kg)				
Size (cm)		21-30		25-30
Low - Moderate Impact	1	2	3	4

CPI Spectra Sock	Size (cm)
NCPI-SS-SM	Small (19-25cm)
NCPI-SS-LG	Large (26-31cm)
NCPI-SSK-SM	Small 6 pack (19-25cm)
NCPI-SSK-LG	Large 6 pack (26-31cm)

Spare Parts Ordering Information

Foot Shell	Size (cm)
NCS-S ■ ● ▲	21-30

■ = "L" or "R" for Left or Right
 ● = Colour (C = Caucasian, T = Tan, B = Brown)
 ▲ = Foot size in cm (21-30)

Shelltread	Size (cm)
NCS-T ■ N ▲	21-30

■ Neutral

Product Specification

Mounting	Sizes	Weight Limit	Build Height
Endo (IP)	21-24cm	100kg	7.8cm
	25-30cm	136kg	8.4cm

Impact Level	Weight of Unit	Warranty	Foot Shell
L M H	744g Size 26 inc. foot shell	2 years Foot shell - 6 months	C T B Caucasian, Tan or Brown Standard width

Technical Information

