

Kinegen Wave



The Kinegen Wave provides users with a number of technical advantages which are necessary for optimal safety and natural movement.

The high-quality knee joint is complemented by the modern design that is seen in the well-known models of the Kinegen series. A polycentric knee joint that has a high stance phase stability due to a posteriorly positioned pivot point, suitable for K2 users. The mechanical shortening in the swing phase enables increased ground clearance whilst the manual lock provides additional safety. Users can move freely and carefree in fresh, salt or chlorinated water thanks to its water-resistant materials.

Waterproof*

*Waterproof note: Following exposure to chlorinated water, rinse immediately with fresh water and dry.



Features and Benefits


- Integrated extension assist spring
- High stance phase safety due to polycentric design
- Ability to easily switch between free knee and locked knee
- 125kg weight limit
- Male pyramid & M36 thread connection options
- Adjustable axis friction for swing phase control
- Waterproof in fresh, chlorinated and salt water*



Adjustable axle friction for swing phase control



Manual lock

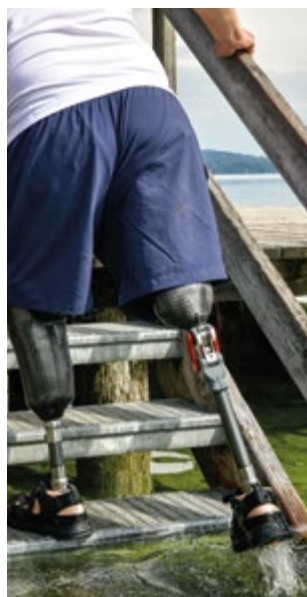

Integrated channels for water drainage



Aluminium/stainless steel

Knees - Streifeneder - Kinegen Wave

Ordering Information

Trial Period - There is no trial period on this knee, however please contact Customer Services or a Steeper Product Manager to enquire.

Impact Level Descriptions

Low: Low impact activities include walking on uneven surfaces.

Kinegen Wave	
Part Number	Description
NMA3A900	KINEGEN Wave with Pyramid adapter
NMA3A900/1	KINEGEN Wave with M36 Thread adapter

Technical Information

	NMA3A900	NMA3A900/1
Proximal Connection	Pyramid adapter	M36 threaded adapter
Activity Level	K2	K2
Impact Level		
Maximum Patient Weight	125kg	125kg
Distal Connection	30mm pylon receiver	30mm pylon receiver
Knee Flexion Angle	145°	145°
Weight	715g	700g
Total Height	165mm	157mm
Material	Aluminium/Stainless Steel	Aluminium/Stainless Steel
Warranty	2 years	2 years

