

The Streifeneder silicone liners offer solutions to meet a wide range of user requirements. Each liner uses robust, first-class materials, alongside Streifeneder's new knit technology, to provide a great balance of comfort and function.

The ClassicSil.basic is a 3mm silicone liner made with firm silicone.



Features and Benefits


- Ideal introductory silicone liner for residual limbs with reasonable soft tissue coverage
- Suitable for low, moderate and high activity levels
- Designed for those looking for improved comfort and suspension
- Lined with antiseptic additives called 'SteriOne', to help keep liners free of bacteria and micro-organisms
- Cost effective 1 year warranty
- Integrated 10cm knitted matrix reduces distal elongation under tensile load
- Easy to don and doff due to resilient polyamide/ elastane cover
- Integrated coloured vertical stripe for easier donning orientation
- Distal silicone-rotation control
- Available with or without distal connections



Antibacterial Technology (SteriOne)




Suitable for Water Use

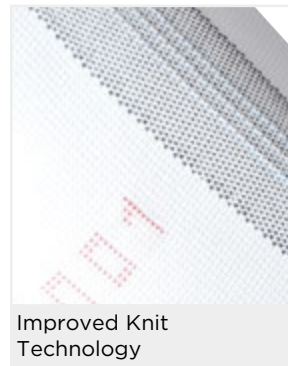

Improved Knit Technology

ComfortSil.basic

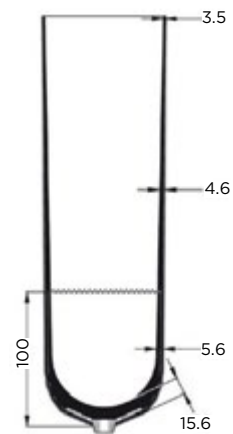

Ordering Information

Ordering Example NMA50L10/18 (with distal connection, size 18)				
Size	Thickness	Outer Matrix	Distal Connection	Part Number
16, 18, 20, 21, 22, 23.5, 25, 26.5, 28, 30, 32, 34, 36, 38, 40, 42, 45	3mm	10cm	with	NMA50L10/Size
16, 18, 20, 21, 22, 23.5, 25, 26.5, 28, 30, 32, 34, 36, 38, 40, 42, 45	3mm	without	without	NMA50L20/Size



**Sizing Note:** Measure the distal circumference approx. 4cm above the distal end and select the liner in the same size or one size smaller.

Suitability Chart

Amputation Type	ClassicSil.basic
Pear-shaped below-knee (Transtibial)	● 3mm basic
Normal, cylindrical below-knee (Transtibial)	●●● 3mm basic

Key: Optimal: ●●● Good: ●● Suitable: ●

Technical Information

Activity Level	K2 to K4
Wall Thickness	3mm
Total Length	Approx 370mm
Material	Firm silicone, medical grade, approx. 40 shore 00
Amputation Level	Transtibial
Distal Connector	Pin Adapter or without
Warranty	1 Year
Colour	Silver Grey

Supplied with a drying stand.