



Uniprox  
SoftSkinAir®

 **steeper**  
by equal.



SoftskinAir® Range - Overview



**SoftSkinAir®  
NS30/NS40/NS50**  
Breathable Silicone Liners

**SoftSkinAir®  
NS30CL/NS40CL/NS50CL**  
Breathable Silicone Liners

**SoftSkinAir®  
2L**  
Breathable Silicone Liners

**SoftSkinAir®  
SPORT**  
Breathable Silicone Liner

**SoftSkinAir®  
Cushion 30** | **SoftSkinAir®  
Cushion 40**  
Breathable Silicone Liners

**SoftSkinAir®  
TF**  
Comfortable & Breathable  
Transfemoral Liners

Material	Silicone (with micropores)				Silicone (with micropores)		Silicone (with micropores)		Silicone (With micropores)		Silicone (with micropores)		Silicone (with micropores)	
Shore Hardness	30/40/50		30/40/50		30/50		40		30		40		40	
Thickness	3mm		3mm		3mm		3mm		3mm		3mm		3mm	
Length	370mm		370mm		370mm		370mm		370mm (23.5 - 28mm) 400mm (28 - 36mm) 430mm (36mm<)		370mm (23.5 - 28mm) 400mm (28 - 36mm) 430mm (36mm<)		350mm	
Amputation Level	Transtibial		Transtibial/Transfemoral		Transtibial		Transtibial		Transtibial/Transfemoral		Transtibial/Transfemoral		Transfemoral	
Distal Connection	Yes		No		Yes		Yes		No		No		Yes	
Mobility Grades	NS30 & NS40 1,2 & 3		NS50 2,3 & 4		NS30CL & NS40CL 1,2 & 3		NS50CL 2,3 & 4		1,2 & 3		1,2 & 3		1,2,3 & 4	
Warranty	6 months		6 months		6 months		6 months		6 months		6 months		6 months	

# SoftSkinAir® NS30/NS30CL


The SoftSkinAir® is a unique, breathable silicone liner featuring a distinctive structure and special micropores to reduce the build up of sweat and draw moisture away from the residual limb; providing a unique solution to patients who previously have been unable to wear a liner due to high sweat levels. It has proven effectiveness in wearing tests and moisture level analyses.

Available in a 30 shore hardness with a distal connector or a smooth distal contour.



## Features and Benefits

- Shore hardness: 30 shore
- Wicks sweat via micropores
- Unique functional and progressive outside texture
- Available with or without distal connector
- Machine-washable (30°C)
- Sizes: 16 – 45
- Breathable and comfortable solution
- Improved grip contact
- Unique outer textile allows for easier socket donning
- Decreased skin irritation
- Reduced temperature inside the liner

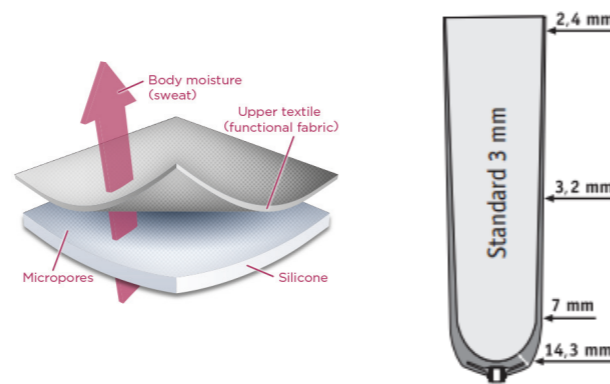

## Ordering Information


All liners are 3mm thick, 370mm long.

**Note:** The circumference is measured 40mm above the distal end of the residual limb

Size	With Distal Connector	Without Distal Connector
	SoftSkinAir NS30	SoftSkinAir NCL30
16	NS30-16-AIR	NS30CL-16-AIR
18	NS30-18-AIR	NS30CL-18-AIR
20	NS30-20-AIR	NS30CL-20-AIR
21	NS30-21-AIR	NS30CL-21-AIR
22	NS30-22-AIR	NS30CL-22-AIR
23.5	NS30-23-AIR	NS30CL-23-AIR
25	NS30-25-AIR	NS30CL-25-AIR
26.5	NS30-26-AIR	NS30CL-26-AIR
28	NS30-28-AIR	NS30CL-28-AIR
30	NS30-30-AIR	NS30CL-30-AIR
32	NS30-32-AIR	NS30CL-32-AIR
34	NS30-34-AIR	NS30CL-34-AIR
36	NS30-36-AIR	NS30CL-36-AIR
38	NS30-38-AIR	NS30CL-38-AIR
40	NS30-40-AIR	NS30CL-40-AIR
42	NS30-42-AIR	NS30CL-42-AIR
45	NS30-45-AIR	NS30CL-45-AIR

## Technical Information

	SoftSkinAir® NS30		SoftSkinAir® NS30CL	
Material	Silicone with micropores		Silicone with micropores	
Shore Hardness	30		30	
Thickness	3mm		3mm	
Length	370mm		370mm	
Amputation Level	Transtibial		Transtibial	
Distal Connector	With distal connector		Without distal connector	
Mobility Grades	NS30 & NS40	NS50	NS30CL & NS40CL	NS50CL
	1,2 & 3	2,3 & 4	1,2 & 3	2,3 & 4
Warranty	6 months		6 months	

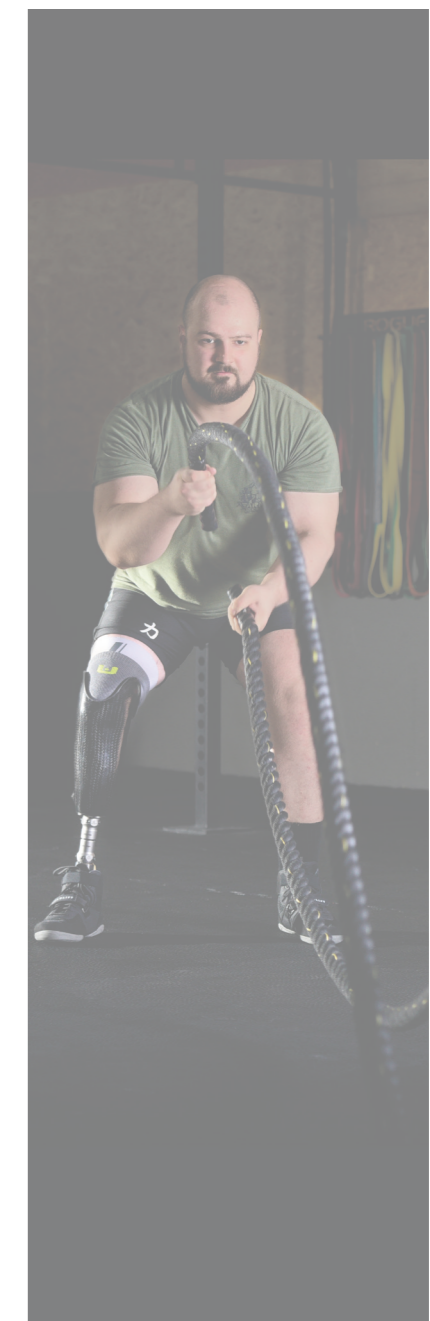

NS30/NS30CL



# SoftSkinAir® NS40/NS40CL


The SoftSkinAir® is a unique, breathable silicone liner featuring a distinctive structure and special micropores to reduce the build up of sweat and draw moisture away from the residual limb; providing a unique solution to patients who previously have been unable to wear a liner due to high sweat levels. It has proven effectiveness in wearing tests and moisture level analyses.

Available in a 40 shore hardness with a distal connector or a smooth distal contour.



## Features and Benefits

- Shore hardness: 40 shore
- Wicks sweat via micropores
- Unique functional and progressive outside texture
- Available with or without distal connector
- Machine-washable (30°C)
- Sizes: 16 – 45
- Breathable and comfortable solution
- Improved grip contact
- Unique outer textile allows for easier socket donning
- Decreased skin irritation
- Reduced temperature inside the liner

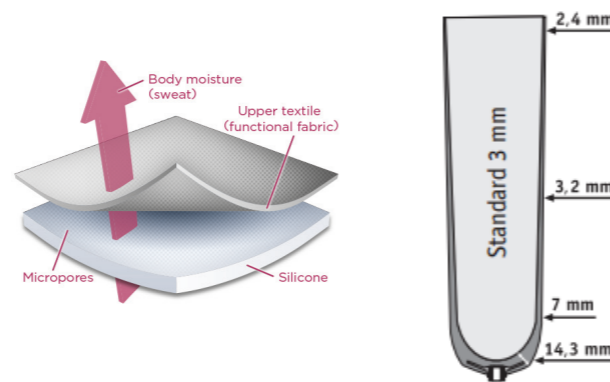

## Ordering Information

All liners are 3mm thick, 370mm long.

**Note:** The circumference is measured 4cm above the distal end of the residual limb

Size	With Distal Connector	Without Distal Connector
	SoftSkinAir NS40	SoftSkinAir NCL40
16	NS40-16-AIR	NS40CL-16-AIR
18	NS40-18-AIR	NS40CL-18-AIR
20	NS40-20-AIR	NS40CL-20-AIR
21	NS40-21-AIR	NS40CL-21-AIR
22	NS40-22-AIR	NS40CL-22-AIR
23.5	NS40-23-AIR	NS40CL-23-AIR
25	NS40-25-AIR	NS40CL-25-AIR
26.5	NS40-26-AIR	NS40CL-26-AIR
28	NS40-28-AIR	NS40CL-28-AIR
30	NS40-30-AIR	NS40CL-30-AIR
32	NS40-32-AIR	NS40CL-32-AIR
34	NS40-34-AIR	NS40CL-34-AIR
36	NS40-36-AIR	NS40CL-36-AIR
38	NS40-38-AIR	NS40CL-38-AIR
40	NS40-40-AIR	NS40CL-40-AIR
42	NS40-42-AIR	NS40CL-42-AIR
45	NS40-45-AIR	NS40CL-45-AIR

## Technical Information

	SoftSkinAir® NS40		SoftSkinAir® NS40CL	
Material	Silicone with micropores		Silicone with micropores	
Shore Hardness	40		40	
Thickness	3mm		3mm	
Length	370mm		370mm	
Amputation Level	Transtibial		Transtibial	
Distal Connector	With distal connector		Without distal connector	
Mobility Grades	NS30 & NS40	NS50	NS30CL & NS40CL	NS50CL
	1,2 & 3	2,3 & 4	1,2 & 3	2,3 & 4
Warranty	6 months		6 months	



NS40/NS40CL



## SoftSkinAir® NS50/NS50CL


The SoftSkinAir® is a unique, breathable silicone liner featuring a distinctive structure and special micropores to reduce the build up of sweat and draw moisture away from the residual limb; providing a unique solution to patients who previously have been unable to wear a liner due to high sweat levels. It has proven effectiveness in wearing tests and moisture level analyses.

Available in a 50 shore hardness, with a distal connector, or with a smooth distal contour.



### Features and Benefits

- Shore hardness: 50 Shore
- Wicks sweat via micropores
- Unique functional and progressive outside texture
- Available with or without distal connector
- Machine-washable (30°C)
- Sizes: 16 – 45
- Breathable and comfortable solution
- Improved grip contact
- Unique outer textile allows for easier socket donning
- Decreased skin irritation
- Reduced temperature inside the liner

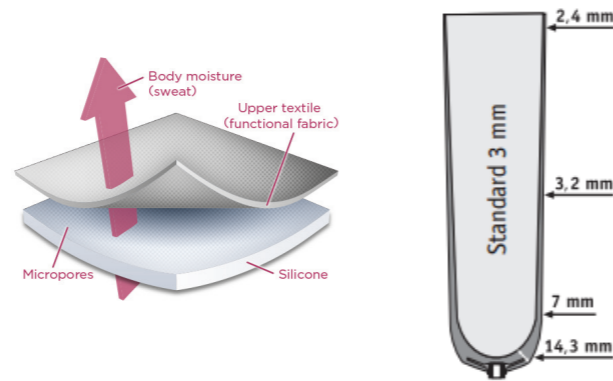

### Ordering Information


All liners are 3mm thick, 370mm long.

**Note:** The circumference is measured 4cm above the distal end of the residual limb

Size	With Distal Connector	Without Distal Connector
	SoftSkin Air NS50	SoftSkin Air NCL50
16	NS50-16-AIR	NS50CL-16-AIR
18	NS50-18-AIR	NS50CL-18-AIR
20	NS50-20-AIR	NS50CL-20-AIR
21	NS50-21-AIR	NS50CL-21-AIR
22	NS50-22-AIR	NS50CL-22-AIR
23.5	NS50-23-AIR	NS50CL-23-AIR
25	NS50-25-AIR	NS50CL-25-AIR
26.5	NS50-26-AIR	NS50CL-26-AIR
28	NS50-28-AIR	NS50CL-28-AIR
30	NS50-30-AIR	NS50CL-30-AIR
32	NS50-32-AIR	NS50CL-32-AIR
34	NS50-34-AIR	NS50CL-34-AIR
36	NS50-36-AIR	NS50CL-36-AIR
38	NS50-38-AIR	NS50CL-38-AIR
40	NS50-40-AIR	NS50CL-40-AIR
42	NS50-42-AIR	NS50CL-42-AIR
45	NS50-45-AIR	NS50CL-45-AIR

### Technical Information

	SoftSkinAir® NS50		SoftSkinAir® NS50CL	
Material	Silicone with micropores		Silicone with micropores	
Shore Hardness	50		50	
Thickness	3mm		3mm	
Length	370mm		370mm	
Amputation Level	Transtibial		Transtibial	
Distal Connector	With distal connector		Without distal connector	
Mobility Grades	NS30 & NS40	NS50	NS30CL & NS40CL	NS50CL
	1,2 & 3	2,3 & 4	1,2 & 3	2,3 & 4
Warranty	6 months		6 months	

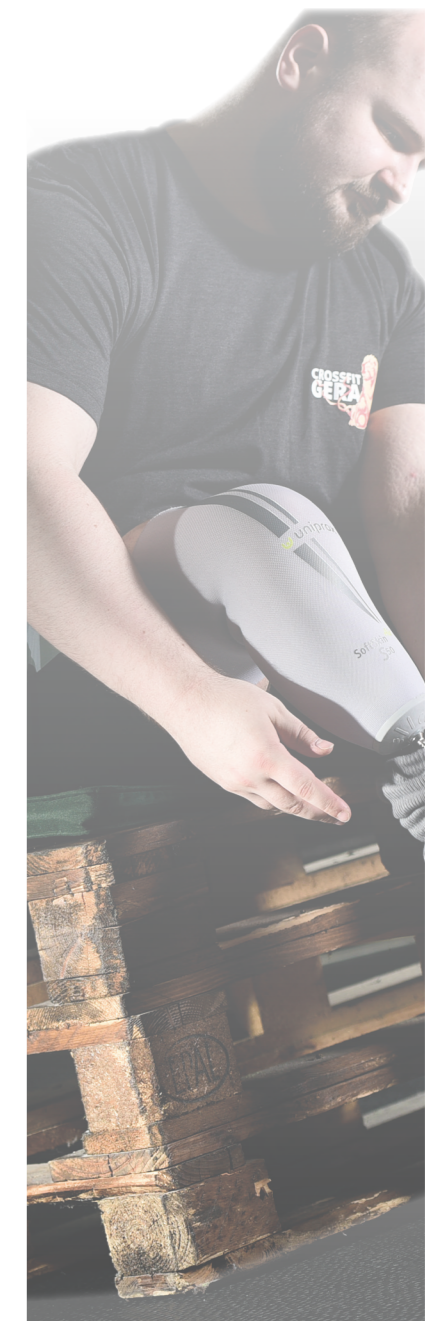

NS50/NS50CL



## SoftSkinAir® 2L

The SoftSkinAir® 2L is a dual layer liner featuring a soft silicone inner for gentle skin contact, with a firmer silicone outer to provide stability.


The liner offers sweat wicking micropores to reduce moisture, resulting in drier skin and lower temperature within the liner for residual limb comfort and reduced skin irritation.

Available in a dual 30 (inner) and 50 (outer) shore hardness, with a distal connector.



### Features and Benefits

- Soft on the inside, stable on the outside (Shore hardness 30/50)
- Soft residual limb padding
- Innovative functional knitted fabric to reduce the longitudinal movement without an internal matrix
- Sizes: 16 - 45
- Unique functional outer fabric
- Machine-washable (30°C)
- Improved grip contact
- Decreased skin irritation
- Breathable and comfortable solution
- Wicks sweat via micropores

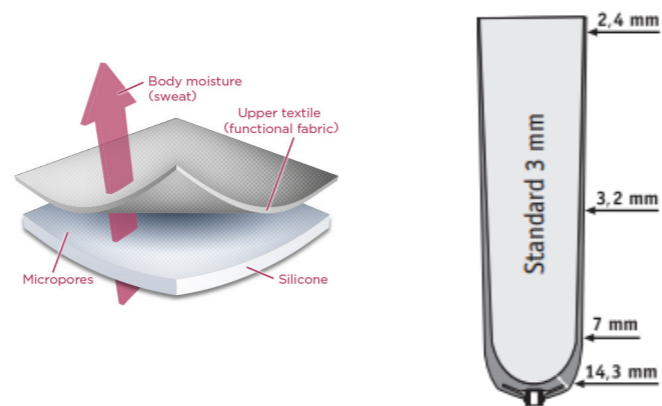

### Ordering Information

All liners are 3mm thick, 370mm long.

**Note:** The circumference is measured 4cm above the distal end of the residual limb

Size	SoftSkinAir® 2L
16	NS-2L-16-AIR
18	NS-2L-18-AIR
20	NS-2L-20-AIR
21	NS-2L-21-AIR
22	NS-2L-22-AIR
23.5	NS-2L-23-AIR
25	NS-2L-25-AIR
26.5	NS-2L-26-AIR
28	NS-2L-28-AIR
30	NS-2L-30-AIR
32	NS-2L-32-AIR
34	NS-2L-34-AIR
36	NS-2L-36-AIR
38	NS-2L-38-AIR
40	NS-2L-40-AIR
42	NS-2L-42-AIR
45	NS-2L-45-AIR

### Technical Information

	SoftSkinAir® 2L
Material	Silicone with micropores
Shore Hardness	30 inner/50 outer
Thickness	3mm
Length	370mm
Amputation Level	Transtibial
Distal Connector	With distal connector
Mobility Grades	1,2 & 3
Warranty	6 months



SoftskinAir® 2L



## SoftSkinAir® Sport


The SoftSkinAir® Sport is a unique, breathable silicone liner offering an additional distal matrix to further improve the reduction of longitudinal movements. This liner features a distinctive structure and special micropores to reduce the build up of sweat and draw moisture away from the residual limb; providing a unique solution to patients who previously have been unable to wear a liner due to high sweat levels.

Available in a 40 shore hardness with a distal connector.



### Features and Benefits

- Shore hardness: 40 shore A (soft)
- Excellent wicking of sweat through micropores
- Unique functional and progressive outside texture
- Available with or without distal connector
- Machine-washable (30°C)
- Breathable and comfortable solution
- Improved grip contact
- Unique outer textile allows for easier socket donning
- Decreased skin irritation
- Reduced temperature inside the liner
- Additional distal matrix for more stability
- Sizes: 20 - 30
- Dual layer with silky inner surface
- Reduction of elongation

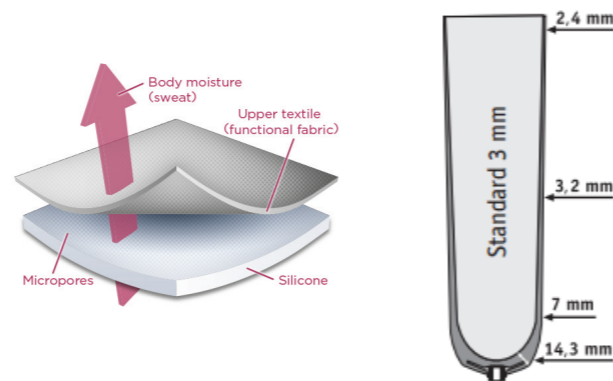

### Ordering Information

All liners are 3mm thick, 370mm long.

**Note:** The circumference is measured 4cm above the distal end of the residual limb

Size	SoftSkinAir® Sport
20	NSPORT-20-AIR
22	NSPORT-22-AIR
23.5	NSPORT-23.5-AIR
25	NSPORT-25-AIR
26.5	NSPORT-26.5-AIR
28	NSPORT-28-AIR
30	NSPORT-30-AIR

### Technical Information

	SoftSkinAir® Sport
Material	Silicone with micropores
Shore Hardness	40
Thickness	3mm
Length	370mm
Amputation Level	Transtibial
Distal Connector	With distal connector
Mobility Grades	2,3 & 4
Warranty	6 months



SoftskinAir® Sport



## SoftSkinAir® Cushion 30 / Cushion 40



The SoftSkinAir® is a unique, breathable silicone liner featuring a distinctive structure and special micropores to reduce the build up of sweat and draw moisture away from the residual limb; providing a unique solution to patients who previously have been unable to wear a liner due to high sweat levels. It has proven effectiveness in wearing tests and moisture level analyses.

Available in 30 and 40 shore hardnesses, with a smooth distal contour.



### Features and Benefits

- Wicks sweat via micropores
- Unique functional and progressive outside texture
- Available with or without distal connector
- Machine-washable (30°C)
- Breathable and comfortable solution
- Improved grip contact
- Unique outer textile allows for easier socket donning
- Decreased skin irritation
- Reduced temperature inside the liner

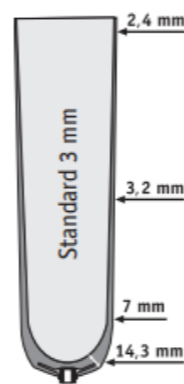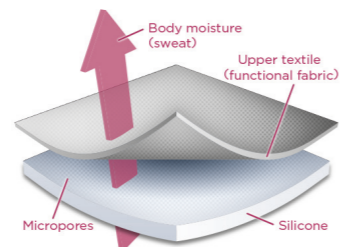

### Ordering Information

All liners are 3mm thick, 370mm long.

**Note:** The circumference is measured 4cm above the distal end of the residual limb

Size	SoftSkinAir® Cushion 30	SoftSkinAir® Cushion 40
23.5	NSC30-23-AIR	NSC40-23-AIR
25	NSC30-25-AIR	NSC40-25-AIR
26.5	NSC30-26-AIR	NSC40-26-AIR
28	NSC30-28-AIR	NSC40-28-AIR
30	NSC30-30-AIR	NSC40-30-AIR
32	NSC30-32-AIR	NSC40-32-AIR
34	NSC30-34-AIR	NSC40-34-AIR
36	NSC30-36-AIR	NSC40-36-AIR
38	NSC30-38-AIR	NSC40-38-AIR
40	NSC30-40-AIR	-
42	NSC30-42-AIR	-
45	NSC30-45-AIR	-

### Technical Information

	SoftSkinAir® Cushion 30	SoftSkinAir® Cushion 40
Material	Silicone with micropores	Silicone with micropores
Shore Hardness	30	40
Thickness	3mm	3mm
Length	370mm (23.5 - 28mm) 400mm (28 - 36mm) 430mm (36mm<)	370mm (23.5 - 28mm) 400mm (28 - 36mm) 430mm (36mm<)
Amputation Level	Transtibial/Transfemoral	Transtibial/Transfemoral
Distal Connector	Smooth contour	Smooth contour
Mobility Grades	1,2 & 3	1,2 & 3
Warranty	6 months	6 months



CUSHION 30 / 40



## SoftSkinAir® TF


The SoftSkinAir® TF silicone liner offers the benefits of the SoftskinAir® range to a wider range of transfemoral users thanks to their increased size range. The TF Air features sweat-wicking micropores to reduce moisture, resulting in drier skin and lower temperature within the liner for residual limb comfort and skin irritation reduction.

With a shore hardness of 40, these liners allow long lasting comfort to the user with high transverse elasticity to adapt to residuum shape.



### Features and Benefits

- Unique transfemoral wicking liner
- Soft residual limb padding
- Perfect sweat evacuation through micropores
- Innovative functional knitted fabric to reduce the longitudinal movement without an internal matrix
- Sizes: 28 - 55
- Unique functional outer fabric
- Machine-washable (30°C)
- Improved grip contact
- Decreased skin irritation
- Breathable and comfortable solution

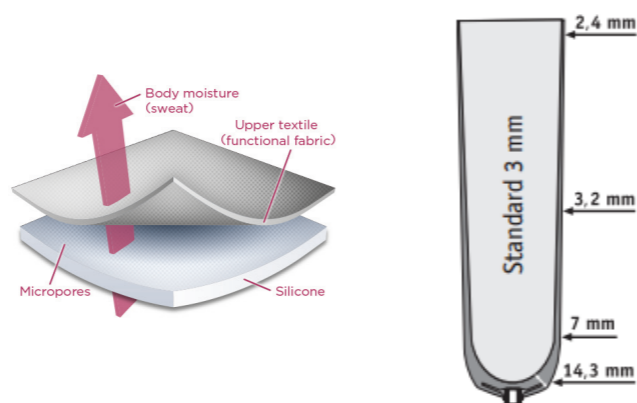

### Ordering Information

All liners are 3mm thick, 370mm long.

**Note:** The circumference is measured 4cm above the distal end of the residual limb

Size	SoftSkinAir® TF
28	NSTF-28-AIR
30	NSTF-30-AIR
32	NSTF-32-AIR
34	NSTF-34-AIR
36	NSTF-36-AIR
38	NSTF-38-AIR
40	NSTF-40-AIR
42	NSTF-42-AIR
45	NSTF-45-AIR
50	NSTF-50-AIR
55	NSTF-55-AIR

### Technical Information

	SoftSkinAir® TF
Material	Silicone with micropores
Shore Hardness	40
Thickness	3mm
Length	350mm
Amputation Level	Transfemoral
Distal Connector	With distal connector
Mobility Grades	1,2,3 & 4
Warranty	6 months

SoftskinAir® TF



Steeper Group  
Unit 3, Stourton Link  
Intermezzo Drive  
Leeds  
LS10 1DF  
United Kingdom



STPPR302  
Tel: +44 (0) 870 240 4133  
Email: [customerservices@steepergroup.com](mailto:customerservices@steepergroup.com)

[www.steepergroup.com](http://www.steepergroup.com)