

Genesee Gel Sleeve

The Genesee Gel Sleeve from College Park can be used as a primary suspension method, or with a suction or vacuum suspension. The sleeve is made of hypo-allergenic material that has been infused with mineral oil to help protect against shear and friction, while also moisturising skin and scar tissue. The sleeve's 3mm thickness provides comfort and knee flexion.

Features and Benefits

- Reinforced interior for puncture-resistance
- Durable fabric-to-gel bond
- Available in 5 sizes
- 2 neutral tones (Grey or Tan)



Reinforced interior




Sleeves- College Park - Genesee Gel Sleeve

Ordering Information

Colour: Grey	
Size	Part Number
1	NGSS-1G
2	NGSS-2G
3	NGSS-3G
4	NGSS-4G
5	NGSS-5G

Colour: Tan	
Size	Part Number
1	NGSS-1T
2	NGSS-2T
3	NGSS-3T
4	NGSS-4T
5	NGSS-5T

Size Selection

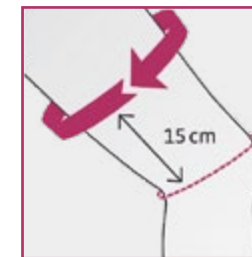

Determining the size
Measure the thigh circumference 15cm above the patella and choose the appropriate size of suspension sleeve.

Technical Information

Proximal Circumference (cm)	Distal Circumference (cm)	Size
22-33	20-30	1
27-38	25-35	2
35-50	30-40	3
40-55	33-48	4
43-60	38-50	5

ContexGel.protect



To avoid premature tearing of the Genesee Gel knee sleeves, we recommend using the ContexGel.protect as a protective cover for sharp-edged prosthetic sockets. This protective cover is made of durable and expandable fabric and has an inner silicone edge on the distal side to prevent slipping.

ContexGel.protect cover				
Size	Side	For Knee Sleeve	Length	Part Number
S	Universal	Size S	19cm	NMA3S50/S
M	Universal	Size M	19cm	NMA3S50/M
L	Universal	Size L	22cm	NMA3S50/L
XL	Universal	Size XL & XXL	22cm	NMA3S50/XL

Protective Inner Sleeve



These sleeve protectors from Össur are made of durable fabric that extends the life of the suspension sleeve and protects it from damage caused by the socket's edge. Silicone pimples on the bottom edge of the sleeve hold it securely in place on the socket.

Protective Sleeve	
Size	Part Number
Small	NM000202
Medium	NM000203
Large	NM000204

Product Specification

Material	Thermoplastic Elastomer (TPE) Gel
Length	38cm
Colours	Grey, Tan
Warranty	3 Months